

Modified 12-Lead ECG

Chest electrodes:

V1 - on the 4th intercostal space at the right sternal border


V2 - on the 4th intercostal space at the left sternal border

V3 - midway between the V2 and V4 electrode positions

V4 - on the 5th intercostal space at the left midclavicular line

V5 - on the left anterior axillary line, horizontal with the V4 electrode position

V6 - on the left midaxillary line, horizontal with the V4 electrode position

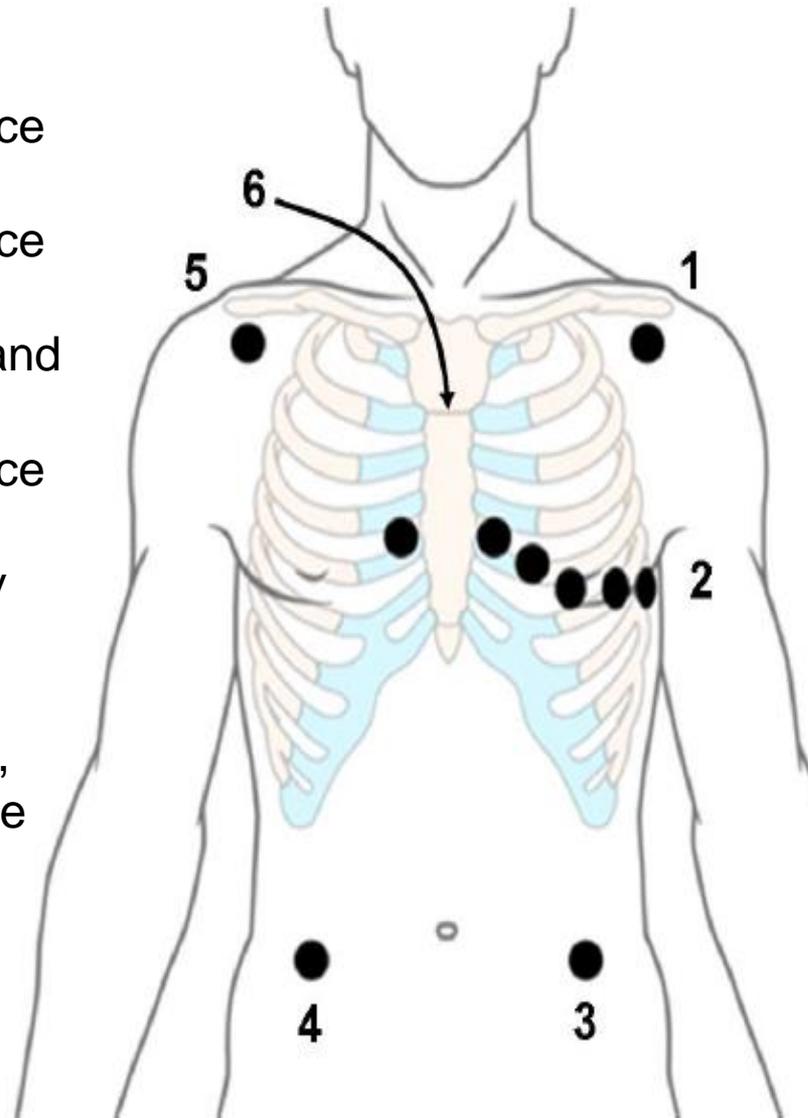

- 1 LA
- 2 V1-V6
- 3 LL
- 4 RL
- 5 RA
- 6 Angle of Louis