


At 0700, the fluid removal resets to "0". You cannot obtain your 0600 to 0700 volume from this screen.

The time represents the fluid removal for the preceding hour. For example, the Periodic volume of 187 ml recorded at 0900 is for the period 0800-0900. The Total volume of 391 ml is the cumulative volume between 0700 and 0900.

The solid white dot to the extreme right indicates this is the current time period.

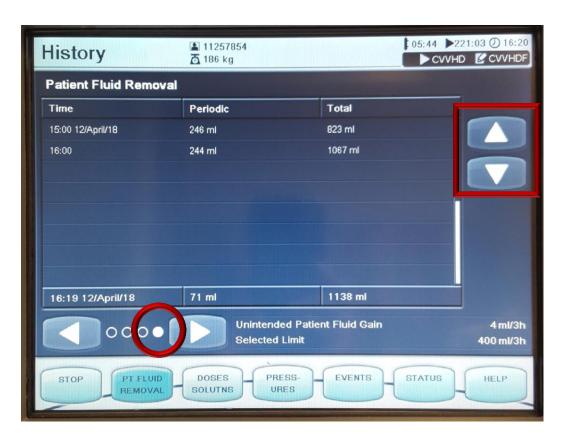

To see other time periods within the current 24 hour period, use the up and down arrows.


The solid white dot to the extreme right indicates that you are scrolling within the current 24 hour period.



To select the period of 0600 to 0700, select the left arrow and go to the previous time period.

23:00 09/April/18 00:00 10/April/18 01:00	246 ml 251 ml	2689 m
00:00 10/April/18	251 ml	2040 0
01:00		2040 11
01:00	240 ml	3180 m
02:00	248 ml	3428 m
03:00	243 ml	3671 m
04:00	248 ml	3919 m
05:00	243 ml	4162 m
06:00	199 ml	4361 m
07:00 10/April/18	195 ml	4556 n

This solid white circle indicates you are on the time period immediately before the current one.

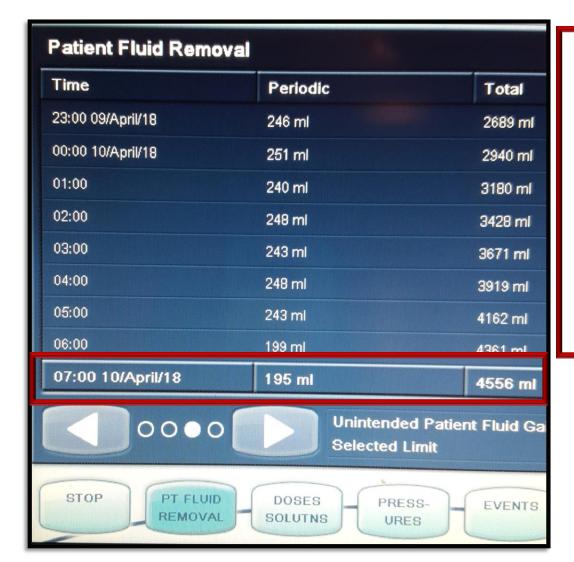

Look at the box immediately below the 0600 time (framed in white).

This is where you will find your 0600 to 0700 fluid removal (Periodic) and the cumulative (Total) volume for the period 0700 to 0700.